



**NEW  
Pediatric  
Sizes!**

**Corbac™  
Adjustable  
Contour  
Back  
Support**



**Forward™  
Cushion**




**Control,  
Comfort,  
Protection**

# It's all about contour...

Whether at work, rest, or play, the Ride Forward Cushion and Corbac surround the pelvis with supportive contours while reducing pressure around at-risk bony prominences. This helps improve postural control without compromise of skin protection.



## Ride Designs® Forward™ Cushion

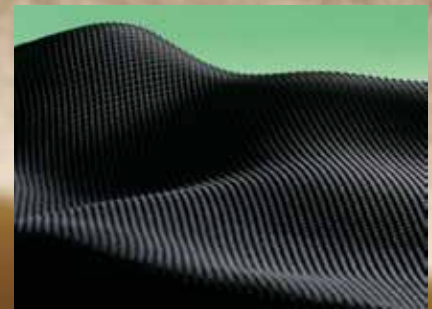

### Control

Ride's patented contoured base helps promote and control improved posture and balance.

### Comfort and protection

Top visco elastic layer provides protection and comfort as it conforms to body shape decreasing peak pressures.

### Even more protection

Even under significant load, the spacer fabric outer cover does not compress. This allows air to move between the sitter and the cushion, helping to keep skin dry and cool.



▲ Spacer fabric in cover.




▲ Inner incontinent resistant cover, shown from the rear.

The cushion contour, spacer fabric outer cover and inner incontinent resistant cover all work together to keep the user and cushion dry. The inner cover really works. There are no seams or zippers in the path of liquid drainage, ensuring the cushion will not absorb moisture. Large accidents or spills simply pour off the back of the cushion.



# Ride Designs® Corbac™ Adjustable Contour Back Support

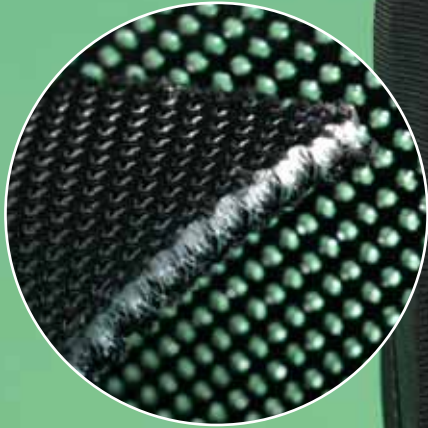

The unique spacer fabric material is used in both the Corbac and Forward Cushion for improved comfort and breathability.



▲ First, dial in the exact posterior contour with the Corbac's adjustment straps.

Then, refine the fit even further with the thoraco-lumbar-pelvic (TLP) supports. TLP supports provide precise adjustment of posterior and lateral contours. ▶

## Lightweight

At two pounds, the Corbac™ weighs little more than standard back upholstery.

## Precision fit

The Corbac is easily adjusted to help improve postural control and alignment. It does not decrease the wheelchair's depth or width. The sitter experiences fitting in, not sitting on, the wheelchair.



## Reduce heat and moisture buildup

The spacer fabric overlay promotes air circulation and wicks moisture from the back surface.

## Versatile

No screws, nuts, or tricky hardware, just a quick and easy mounting system that takes seconds to assemble onto the wheelchair.

## A cinch to clean

The overlay detaches easily and can be machine washed.

# Specs



**NEW!**  
**P**  
 12 x 12"  
 Forward Cushion,  
 contoured & sized for  
 pediatric application.  
 12" Corbac also  
 available.

## Ride Forward Cushion Sizes

Size (width x depth)	Part Number	Fits Wheelchair	
		Widths	Depths
12 x 12" <b>P</b>	FC-1212	11/12"	12/13"
14 x 14"	FC-1414	13/14"	14/15"
14 x 16"	FC-1416	13/14"	16/17"
16 x 16"	FC-1616	15/16"	16/17"
16 x 18"	FC-1618	15/16"	18/19"
16 x 20"	FC-1620	15/16"	20/21"
18 x 16"	FC-1816	17/18"	16/17"
18 x 18"	FC-1818	17/18"	18/19"
18 x 20"	FC-1820	17/18"	20/21"
20 x 16"	FC-2016	19/20"	16/17"
20 x 18"	FC-2018	19/20"	18/19"
20 x 20"	FC-2020	19/20"	20/21"

## Ride CAM® Wedge Kit for Forward Cushion (FCC-WK)

Available for every cushion size.

Use CAM Wedges **symmetrically** to promote more upright posture and decrease ischial pressure.



Use CAM Wedges **asymmetrically** to help level a flexible pelvic obliquity.



## Ride Corbac Sizes

Wheelchair Width	Corbac Height	Part Number
12"	12 - 14"	CB-1200
14"	12 - 15"	CB-1400
15"	13 - 16"	CB-1500
16"	14 - 17"	CB-1600
16"	17 - 20"	CB-1600T
17"	14 - 17"	CB-1700
18"	14 - 17"	CB-1800
18"	17 - 20"	CB-1800T
19 - 20"	14 - 17"	CB-1920
19 - 20"	17 - 20"	CB-1920T

▶ Note: For people at higher risk for skin breakdown, and/or more challenging postural control needs, consider the **Ride Custom Cushion and Back**, the ultimate in postural support and skin care.



**Ride Designs®**  
 a branch of Aspen Seating, LLC  
 4211-G South Natches Court  
 Sheridan, Colorado 80110

Protected by one or more of the following US  
 Patents: 7,373,678; 7,216,388; 7,220,376;  
 6,990,744; 7,395,566; 7,140,057.  
 © 2010, Ride Designs. 7.10 Rev. C.  
 BRO005

phone: 303.781.1633  
 toll-free: 866.781.1633  
[www.ridedesigns.com](http://www.ridedesigns.com)



**PARAGON**  
 mobility

[www.paragonmobility.com.au](http://www.paragonmobility.com.au)  
[info@paragonmobility.com.au](mailto:info@paragonmobility.com.au)